

Respyr Metabolism 1.0 Quick Start Guide


BREATH-BASED METABOLIC WELLNESS DEVICE

Respyr is a trademark of Humors Tech Pvt. Ltd.



Humors Tech Private Limited
2nd Floor, No. 14/7, 3rd Cross, Sunleo Power
Networks, Old Airport Road, Konena Agrahara,
Bengaluru, Karnataka, 560017, India
Contact: +91-82963 80628
Website: www.respyr.us

Version: UM/RESPYRMETABOLISM/01
Date: 23-Feb-2026

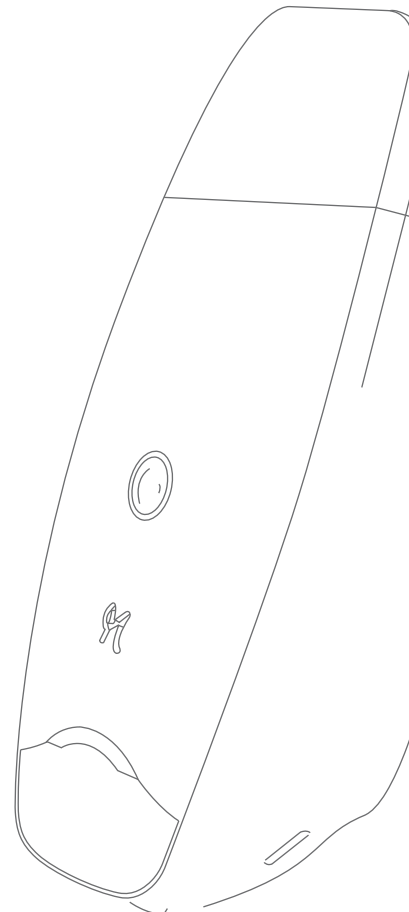

Breath-Based Metabolic Wellness Device

Respyr™ is a handheld breath analysis device designed to measure specific breath biomarkers and provide metabolic wellness insights through the Respyr mobile application.

Please read this manual carefully before using the device and keep it for future reference.

What’s in the Box:

Your Respyr™ package contains the following items illustrated in figure 1:

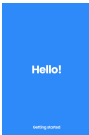
Item No.	Visual	Component	Description
1.		Respyr Device	Handheld breath analysis device used to perform breath measurements.
2.		Charging Dock with Cable	Used to safely charge the Respyr device using a compatible USB power adapter.
3.		Travel Pouch	Used to safely store and carry the Respyr device.
4.		Respyr App Start Guide	Card containing QR codes to download the Respyr mobile application and access the full Instructions for Use (IFU).

fig 1: Respyr package content

Steps:

Steps	Description
Step-1	Download the Respyr app. Scan the QR code on the App Card and install the app on your phone.
Step-2	First charge is critical. Switch ON the device, then place it in the charging dock. Keep the device switched ON throughout the entire charge (approx. 2-3 hrs). <i>After storage of 1 month or more, always perform this full charge before use</i>
Step-3	Associated with digestive activity and gut-related metabolic processes.

Respyr™ must be charged using the provided charging dock and cable connected to a compatible USB power adapter.

First-Time Use Recommendation

Fully charge the Respyr™ device before first use, or after one month or more of inactivity, to ensure proper operation and battery performance.

 **Important:** Switch the device ON before placing it on the charging dock, and keep it ON until charging is complete. This ensures the sensors fully activate and deliver accurate readings.

Steps to Charge Respyr Device:

Steps	Instruction
*	First-Time Use and After Long Inactivity: Switch ON the device before placing it in the charging dock. Keep the device switched ON throughout the entire charging process until the dock indicator light turns OFF.
1.	Connect the charging dock cable to a USB power adapter (5 V DC, minimum 1 A).
2.	Plug the USB power adapter into a power outlet.
3.	Place the Respyr device into the charging dock properly.
4..	Ensure the device is seated correctly in the dock.
5..	Verify that the Red Dot indicator on the dock turns ON, indicating charging has started.
6..	Allow the device to charge until the dock indicator light turns OFF, indicating full charge.

Steps to Perform Breath Test:

Steps	Instruction
1.	Press the power button (short press) to switch ON the Respyr device.
2.	Verify that the Blue LED appears, indicating the device is powered ON.
3.	Open the Respyr mobile application.
4..	Connect the Respyr device to the mobile application via Bluetooth.
5..	Follow the breathing instructions displayed in the mobile application.
6..	Exhale into the Respyr device as instructed by the application.
7.	Keep the device steady during the measurement.
8.	Wait for the measurement to complete and view the results in the mobile application.

Use indoors only, in a clean and well-ventilated environment. Respyr™ is a wellness device — not for medical or clinical use.

Respyr™ uses LED lights to indicate device status and operating conditions. The LED color and behavior help you understand the current state of the device.

Device LED Indicators:

LED Color	Light Behavior	Device Status	Description
Blue	Solid	Device Power ON	Appears when the power button is pressed and the device is switched ON. The device is starting and ready for operation.
Green	Solid	Connected to Respyr App	Indicates that the device is successfully connected to the Respyr mobile application via Bluetooth.
Blue	Solid	Measurement in Progress	Indicates that the device is performing breath measurement.
Red	Blinking	Low Battery	Indicates that the battery level is low. Please charge the device soon.
Red	Slow glowing ON and OFF	Critical Battery Level	Indicates that the battery level is too low and the device cannot operate. Immediate charging is required.
Red	Solid	Reading Failure	Indicates that the breath reading was not successful due to improper inhale, exhale, or breath hold. Please follow the instructions in the Respyr mobile application and repeat the test.
Red	Solid	Device Fault or Error	Indicates a device fault or internal error. Restart the device. If the issue continues, contact Humors Tech Pvt. Ltd. support.

NO.	Description
1.	Clean the breath inlet gently with a dry cotton swab after each use. Do not use liquids or solvents.
2.	Store the device in the provided travel pouch when not in use.
3.	Keep away from moisture, direct sunlight, and extreme temperatures.
4.	Do not drop, crush, or expose the device to water.
5.	Battery is built-in and non-replaceable. Contact support if battery runtime drops significantly.