

Respyr Metabolism 1.0 User's Manual


BREATH-BASED METABOLIC WELLNESS DEVICE

Respyr is a trademark of Humors Tech Pvt. Ltd.



Humors Tech Private Limited
2nd Floor, No. 14/7, 3rd Cross, Sunleo Power
Networks, Old Airport Road, Konena Agrahara,
Bengaluru, Karnataka, 560017, India
Contact: +91-82963 80628
Website: www.respyr.us

Version: UM/RESPYRMETABOLISM/01
Date: 23-Feb-2026

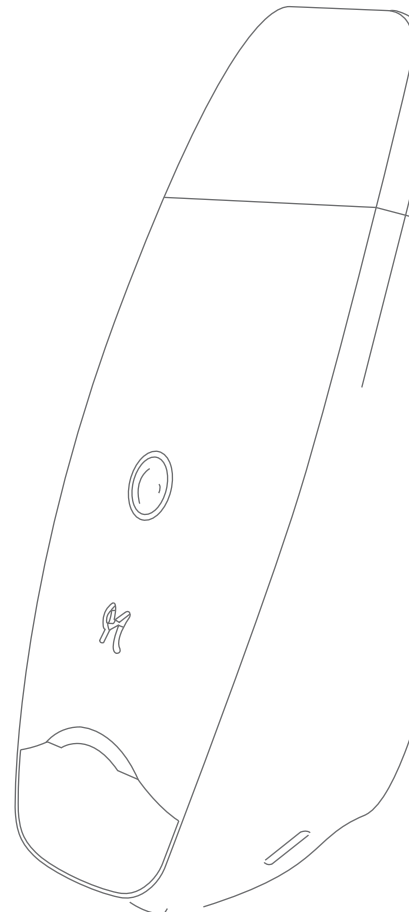

Sl No.	Contents	Page No.
1.	Respyr™ 1.0	01
2.	Contents	02
3.	What's in the Box	03
4.	Device Description	04
5.	LED Indicator Guide	06
6.	Battery and Charging Information	07
7.	Battery Status Indicator	09
8.	Charging Procedure	11
9.	How to Perform a Breath Test	12
10.	Safety Instructions and Warnings	14
11.	Cleaning and Maintenance Instructions	16
12.	Storage and Transport Instructions	17
13.	Manufacturer Information and Support	19
14.	Troubleshooting, Warranty, Return and Refund Policy	21
15.	Copyright Notice and Trademarks	24
16.	Manufacturer Responsibility and Disclaimer	24

Breath-Based Metabolic Wellness Device

Respyr™ is a handheld breath analysis device designed to measure specific breath biomarkers and provide metabolic wellness insights through the Respyr mobile application.

Please read this manual carefully before using the device and keep it for future reference.

What’s in the Box:

Your Respyr™ package contains the following items illustrated in figure 1:

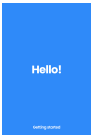
Item No.	Visual	Component	Description
1.		Respyr Device	Handheld breath analysis device used to perform breath measurements.
2.		Charging Dock with Cable	Used to safely charge the Respyr device using a compatible USB power adapter.
3.		Travel Pouch	Used to safely store and carry the Respyr device.
4.		Respyr App Start Guide	Card containing QR codes to download the Respyr mobile application and access the full Instructions for Use (IFU).

fig 1: Respyr package content

Respyr™ is a handheld electronic device designed to analyze exhaled breath and provide metabolic wellness insights through the Respyr mobile application. Respyr measures specific volatile organic compounds (VOCs) present in exhaled breath. These breath biomarkers are associated with metabolic and physiological processes in the body.

The measured breath biomarkers include:

Biomarker	Description
Acetone	Associated with fat metabolism and energy utilization.
Ethanol	Associated with metabolic and biochemical processes in the body.
Hydrogen	Associated with digestive activity and gut-related metabolic processes.

Respyr measures the concentration patterns of these breath biomarkers using integrated sensors. The device transmits the breath data to the Respyr mobile application via Bluetooth, where the data is analyzed and presented as wellness insights.

Intended Use and User Population:

Respyr™ is intended for use by healthy adults seeking to monitor metabolic wellness through breath analysis. The device is designed for personal, non-clinical, indoor use only.

Intended User: Healthy adults (18 years and above) for personal wellness monitoring.

Not intended for: Children, individuals with respiratory conditions, or for clinical diagnostic purposes. This device is not a medical device and does not provide medical diagnoses.

Consult your doctor before using this device if you have any existing medical condition or health concern.

Device Components:



fig 2: Respyr Device

Respyr™ uses LED lights to indicate device status and operating conditions. The LED color and behavior help you understand the current state of the device.

Device LED Indicators:

LED Color	Light Behavior	Device Status	Description
Blue	Solid	Device Power ON	Appears when the power button is pressed and the device is switched ON. The device is starting and ready for operation.
Green	Solid	Connected to Respyr App	Indicates that the device is successfully connected to the Respyr mobile application via Bluetooth.
Blue	Solid	Measurement in Progress	Indicates that the device is performing breath measurement.
Red	Blinking	Low Battery	Indicates that the battery level is low. Please charge the device soon.
Red	Slow glowing ON and OFF	Critical Battery Level	Indicates that the battery level is too low and the device cannot operate. Immediate charging is required.
Red	Solid	Reading Failure	Indicates that the breath reading was not successful due to improper inhale, exhale, or breath hold. Please follow the instructions in the Respyr mobile application and repeat the test.
Red	Solid	Device Fault or Error	Indicates a device fault or internal error. Restart the device. If the issue continues, contact Humors Tech Pvt. Ltd. support.

Dock Light	Light Behavior	Charging Status	Description
Red Dot	Solid	Charging in Progress	Indicates that the Respyr device is currently charging.
No Light	OFF	Fully Charged	Indicates that the device is fully charged.

Battery and Charging Information

Respyr™ contains a built-in rechargeable Nickel-Metal Hydride (NiMH) battery designed to support normal device operation.

Battery Specifications:

Parameter	Specification
Battery Type	Rechargeable Nickel-Metal Hydride (NiMH)
Battery Configuration	2-cell battery pack
Battery Capacity	750 mAh
Nominal Voltage	2.4 V DC
Charging Method	Charging dock with attached cable
Charging Input	5 V DC via charging dock
Charging Time	Couple of hours
Battery Replacement	Battery is not user-replaceable

Note : The power adapter is not supplied with the Respyr device.

Parameter	Requirement
Power Source	USB power adapter
Output Voltage	5 V DC
Output Current	1 A
Safety Certification	Adapter should preferably be UL certified or equivalent safety certified
Charging Dock	Use only the charging dock and cable supplied by Humors Tech Pvt. Ltd.

Charging Safety Instructions:

- Use only the provided charging dock and cable supplied by Humors Tech Pvt. Ltd.
- Use a USB power adapter that is preferably UL certified or safety certified.
- Do not use the device while it is placed in the charging dock.
- Always remove the device from the charging dock before performing a breath measurement.
- Do not use damaged charging adapters or cables.
- Charge the device indoors in a dry and safe environment.
- When using the device for the first time, or after a long period of inactivity (e.g. one month or more), fully charge the device before use. Switch ON the device before placing it in the charging dock and keep it switched ON throughout the entire charging process.


Respyr™ indicates battery status through both the device LED indicator and the Respyr mobile application.

Device Battery Status Indicators:

LED Color	Light Behavior	Battery Status	Description	User Action
Red	Blinking	Low Battery	Indicates that the battery level is becoming low. The device can still operate, but charging is recommended.	Charge the device soon.
Red	Gradually glowing ON and OFF	Critical Battery Level	Indicates that the battery level is extremely low. The device will not allow breath measurement.	Charge the device immediately.
No Red Light	Normal Operation	Sufficient Battery	Indicates that the battery level is sufficient for normal operation.	No action required.

Mobile Application Battery Status Indicator:

App Indicator	Battery Status	Description	User Action
Popup Message: "Low Battery and Charging Required"	Critical Battery Level	When the battery becomes extremely low, the Respyr mobile application displays a popup message stating "Low Battery and Charging Required." The application will not allow the user to perform a breath measurement until the device is charged.	Charge the device using the charging dock before performing a reading.

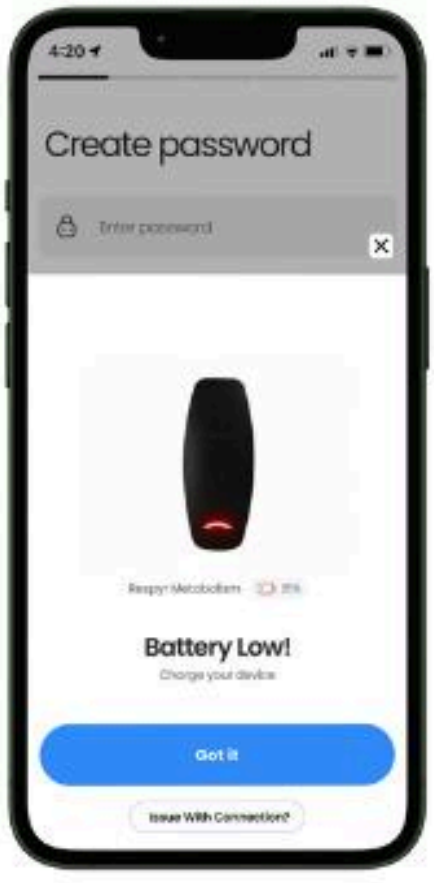

Respyr™ must be charged using the provided charging dock and cable connected to a compatible USB power adapter.

First-Time Use Recommendation

Fully charge the Respyr™ device before first use, or after one month or more of inactivity, to ensure proper operation and battery performance.

 **Important:** Switch the device ON before placing it on the charging dock, and keep it ON until charging is complete. This ensures the sensors fully activate and deliver accurate readings.

Steps to Charge Respyr Device:

Steps	Instruction
*	First-Time Use and After Long Inactivity: Switch ON the device before placing it in the charging dock. Keep the device switched ON throughout the entire charging process until the dock indicator light turns OFF.
1.	Connect the charging dock cable to a USB power adapter (5 V DC, minimum 1 A).
2.	Plug the USB power adapter into a power outlet.
3.	Place the Respyr device into the charging dock properly.
4..	Ensure the device is seated correctly in the dock.
5..	Verify that the Red Dot indicator on the dock turns ON, indicating charging has started.
6..	Allow the device to charge until the dock indicator light turns OFF, indicating full charge.

Respyr™ performs breath measurements through guided instructions provided in the Respyr mobile application. Follow the steps below to perform a correct breath test.

Before Performing the Test:

Instruction	Description
Indoor Ventilated Environment	Perform the test in a clean, indoor, and well-ventilated environment.
Charge Device	Ensure the Respyr device has sufficient battery before performing the test.
Install and Open App	Install and open the Respyr mobile application using the QR code provided in the Respyr App Start Guide.
Enable Bluetooth	Ensure Bluetooth is enabled on your mobile device.

Steps to Perform Breath Test:

Steps	Instruction
1.	Press the power button (short press) to switch ON the Respyr device.
2.	Verify that the Blue LED appears, indicating the device is powered ON.
3.	Open the Respyr mobile application.
4..	Connect the Respyr device to the mobile application via Bluetooth.
5..	Follow the breathing instructions displayed in the mobile application.
6..	Exhale into the Respyr device as instructed by the application.
7.	Keep the device steady during the measurement.
8.	Wait for the measurement to complete and view the results in the mobile application.

Automatic Power OFF:

Condition	Description
After Measurement	The Respyr device automatically switches OFF after completing the breath measurement.

Manual Reset and Power Control:

Action	Instruction
Power ON	Press the power button once (short press). The Blue LED will appear indicating the device is ON.
Hard Reset	Press and hold the power button for 2 seconds to reset the device.

Important Instructions During Measurement:

Instruction	Description
Follow App Guidance	Always follow the breathing instructions provided in the mobile application.
Keep Device Steady	Do not move the device during measurement.
Do Not Block Air Vent	Ensure the air exhaust vent is not blocked.
Repeat if Needed	If reading fails, follow the instructions in the app and repeat the test.

Read and follow all safety instructions before using Respyr™ to ensure safe and proper operation.

General Safety Instructions:

Instruction	Description
Intended Use	Use Respyr only as described in this manual and the Respyr mobile application.
Indoor Use Only	Use the device only in indoor environments.
Keep Device Dry	Do not expose the device to water, liquids, or excessive moisture.
Avoid Physical Damage	Do not drop, crush, or apply excessive force to the device.
Do Not Disassemble	Do not open or disassemble the device. There are no user-serviceable parts inside.

Battery and Charging Safety:

Instruction	Description
Do Not Use While Charging	Do not use the device while it is placed in the charging dock.
Avoid Heat Exposure	Do not expose the device to high temperatures or direct sunlight.
Do Not Damage Battery	Do not puncture, crush, or expose the device to fire.
Non-Replaceable Battery	The Respyr™ device contains a built-in rechargeable battery that is not user-replaceable.
Battery Service	In the event that a significant reduction in battery runtime is observed over time, please contact Humors Tech Pvt. Ltd. for device service or replacement.
Do Not Attempt Replacement	Do not attempt to open the device or replace the battery yourself. This may damage the device and void the warranty.

Note : The battery inside Respyr™ is not user-replaceable. If you notice a significant decrease in battery runtime, please contact Humors Tech Pvt. Ltd. for assistance or replacement of the device.

Device Handling and Storage:

Instruction	Description
Proper Storage	Store the device in a clean and dry place.
Avoid Extreme Conditions	Do not expose the device to extreme temperature or humidity.
Use Travel Pouch	Use the provided travel pouch to protect the device during transport.

Service and Support:

Instruction	Description
No User Serviceable Parts	Do not attempt to repair the device yourself.
Contact Manufacturer	Contact Humors Tech Pvt. Ltd. support for any service or repair needs.

Regular cleaning and proper maintenance help ensure safe operation and reliable performance of the Respyr™ device.

Cleaning Instructions:

Instruction	Description
External Cleaning	Clean the outer surface of the Respyr device using a soft, dry cloth.
Avoid Liquid Entry	Do not allow liquids to enter the breath inlet, air vent, or charging contacts.
Do Not Immerse	Do not immerse the device in water or any liquid.
Do Not Spray	Do not spray cleaning liquid onto the device.

Maintenance Instructions:

Instruction	Description
No User Serviceable Parts	The Respyr device does not contain any user-serviceable parts.
Do Not Open Device	Do not open or attempt to disassemble the device.
Protect from Damage	Avoid dropping or exposing the device to physical impact.
Proper Storage	Store the device in a clean, dry environment when not in use.
Use Travel Pouch	Use the provided travel pouch to protect the device during transport and storage.

Proper storage and transport of the Respyr™ device help maintain device performance and extend product life.

Storage Conditions:

Storage Temperature	_____	0°C to 50°C
Storage Humidity	_____	Up to 90% non-condensing
Storage Environment	_____	Clean, dry, indoor environment
Protection	_____	Store the device in the provided travel pouch when not in use

Storage Instructions:

Instruction	Description
Clean Before Storage	Ensure the device is clean and dry before storing.
Avoid Moisture	Do not store the device in humid or wet environments.
Avoid Direct Sunlight	Do not expose the device to direct sunlight for extended periods.
Avoid Extreme Temperatures	Do not store the device in extremely hot or cold locations.
Safe Placement	Store the device in a safe location to prevent accidental damage.

Instruction	Description
Use Travel Pouch	Always use the provided travel pouch during transport.
Avoid Physical Impact	Protect the device from drops, vibration, or mechanical shock.
Avoid Harsh Conditions	Do not expose the device to excessive heat, cold, or moisture during transport.
Power OFF During Transport	Ensure the device is switched OFF before transport.

Respyr™ is manufactured and supported by Humors Tech Pvt. Ltd.

Manufacturer Details:

Parameter	Information
Manufacturer Name	Humors Tech Pvt. Ltd.
Product Name	Respyr Metabolism
Model	Respyr 1.0
Country of Manufacture	India
Intended Use	Breath-based metabolic wellness device for indoor use

Customer Support:

For technical support, service, or product-related inquiries, please contact:

Contact Type	Details
Website	www.respyr.ai
Email	connect@respyr.in
Support Availability	Contact via email for assistance, service, or replacement requests

Instruction	Description
Authorized Service Only	Device service and repair must be performed only by Humors Tech Pvt. Ltd. or authorized personnel.
Do Not Attempt Repair	Do not attempt to open, repair, or modify the device.
Battery Service	If battery performance significantly decreases, contact Humors Tech Pvt. Ltd. for service or replacement.

Notice	Description
Proper Use	Use the device only as described in this manual.
Safety Compliance	Follow all safety instructions to ensure safe and proper operation.

For troubleshooting assistance, warranty information, return policy, and refund policy, please refer to the official Respyr support pages available online.

These pages provide the latest and most accurate information regarding device usage, support, service, and applicable policies.

You can access the information using the links below:

Troubleshooting Guide:
[connect@respyr.in]

Warranty Policy:
[connect@respyr.in]

Return and Refund Policy:
[connect@respyr.in]

Alternatively, you may scan the QR code provided in the Respyr App Start Guide to access these resources.

For additional support, please contact Humors Tech Pvt. Ltd.:

Website: www.respyr.ai

Email: connect@respyr.in

Device Symbols

The following symbols may appear on the Respyr™ device, packaging, or accompanying documentation. Each symbol is explained in the table below.

Symbol	Meaning
 Serial Number	Identifies the product's unique serial number
 Batch Code	Indicates the production batch
 Manufacturer	Indicates the legal manufacturer
 IVD	In Vitro Diagnostic Device
 Self-testing Device	Intended for self-use
 Date of Manufacture	Indicates the date the product was manufactured
 Temperature Limit	Indicates acceptable storage temperature range
 Biological Risk	Warns of biological contamination risks
 Caution	General caution or warning
 Consult Instructions	Refer to printed user manual before use
 Consult e-IFU	Refer to electronic instructions for use (if applicable)
 ISO 13485:2016	Certified quality management for medical devices
 Disposal	Dispose of device responsibly per local regulations

Symbol	Meaning
 Non-sterile	To indicate that the device that is normally provided sterile in the same or similar packaging has not been sterilized.
 Keep away from rain	To indicate that the transport package shall be kept away from rain and in dry conditions.
 Keep out of reach of children	To indicate packages containing devices and equipment that shall be kept out of reach of children, or to identify devices and equipment that shall be kept out of reach of children.

Safety Warnings and Cautions:

Read all warnings and cautions before using the Respyr™ device. Failure to observe these notices may result in injury or device damage.

WARNING — Indicates a potentially hazardous situation which, if not avoided, may result in personal injury or device damage.

CAUTION — Indicates important information that must be followed to ensure correct and safe use of the device.

NOTE — Indicates supplementary information that may be helpful to the user.

© 2026 Humors Tech Pvt. Ltd. All rights reserved.

Respyr™ and the Respyr logo are trademarks of Humors Tech Pvt. Ltd.

No part of this document may be reproduced, distributed, transmitted, or stored in any form or by any means, electronic, mechanical, photocopying, recording, or otherwise, without prior written permission from Humors Tech Pvt. Ltd.

Humors Tech Pvt. Ltd. reserves the right to make changes to the product and this document without prior notice.

Manufacturer Responsibility and Disclaimer

Humors Tech Pvt. Ltd. is responsible for the safety, performance, and reliability of the Respyr™ device only when:

- The device is used in accordance with the instructions provided in this manual.
- The device is not modified, altered, or repaired by unauthorized personnel.
- Only approved accessories and charging components are used.

Humors Tech Pvt. Ltd. shall not be liable for any damage, malfunction, or injury resulting from:

- Improper use of the device
- Use outside the specified operating conditions
- Unauthorized modification or repair
- Use of unapproved accessories or charging equipment

Respyr™ is intended for use by healthy individuals wellness monitoring purposes only and is not intended to diagnose, treat, cure, or prevent any disease.

The information and results provided by Respyr™ and the Respyr mobile application are intended for informational and wellness purposes only.

Users should always consult a qualified healthcare professional or doctor for medical advice, diagnosis, or treatment regarding any health condition or concern.

Do not use Respyr™ as a substitute for professional medical advice or medical examination.

If you have any medical condition or health concern, consult your doctor before making any health-related decisions.